

*This medical device is intended for patient positioning.
This is not a restraint system.*

*Before use, it is important to check the status of the device: stitching, state of the buckles, slings,
fabric. Any damage will either be repaired by NAUSICAA Medical, or will result in a disposal.*

INTEGRAL BIRDY BELT

FOR ALL ARMCHAIRS



THIS MEDICAL DEVICE MUST BE
IMPLEMENTED BY A TRAINED
HEALTH PROFESSIONAL




**PREVENT FORWARD
SLIDING AND MAINTAIN
THE CHEST IN UPRIGHT
POSITION FOR ALL CHAIRS**



Maximum useful scope: + 220 cm
Buckle: 50 mm

It can be set up with the patient seated and avoids any pelvic constriction.
Its shoulder heights are separately adjustable to help position and compensate for a bad
position.

- Available in 4 sizes S, M, L, XL
- Waist (in cm)

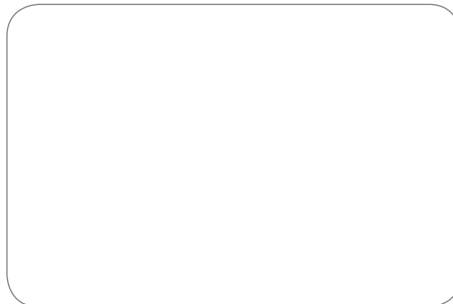

CBRINHP-PF

Important Recommendations :

In order to get a maximum efficiency of these sort of product, it is necessary :

- to choose the appropriate size for the patient
- to get the best possible adjustment to the patient

These products must not be in direct contact with a wounded skin.

INTEGRAL BIRDY BELT

HOSPITAL RANGE
Combination of polyurethan foamand
on both faces
Jersey/PU coating
100% polyester bias
100% polyester strap
Nylon buckle
Thermo compressed foam

INSTRUCTIONS FOR USE

1. Slip the head through the shoulders harness.
2. Place the Birdy belt in the patient's back at the lumbar and adjust the small loops on each hip.
3. Pass the end of the belt under the slightly raised thigh.
4. Do the same for the second thigh.
5. Pass each buckle in the corresponding loop.
6. Set up the abdominal belt with the velcro strap.
7. Place the double buckle in the back strap.
8. Close the buckle at the back of the chair.
9. At the back, set up the two rear buckles of the harness.



v 02.03/2019



THIS MEDICAL DEVICE MUST BE SET UP BY A PROFESSIONAL
OR THE USER MUST BE EDUCATED TO SET IT UP.