

# NAUSICAA WEIGHING SYSTEM



## Easily weigh patients during transfers

All our patient lifts can be equipped with a weighing system which is designed to easily determine the weight of patients (up to 300 kg) during transfers by a patient lifter.


## Good to know

The device's maximum capacity on which it will be assembled prevails on the weighing system.

It can only be mounted on a NAUSICAA Medical patient lift and is delivered with fixations to reduce the loss of useful height.

## Caution

To maintain a correct useful height, we do not recommend using it with compact devices such as the BLUE FLYER COMPACT, the NAUSI-FLY COMPACT and the FLYER COMPACT.

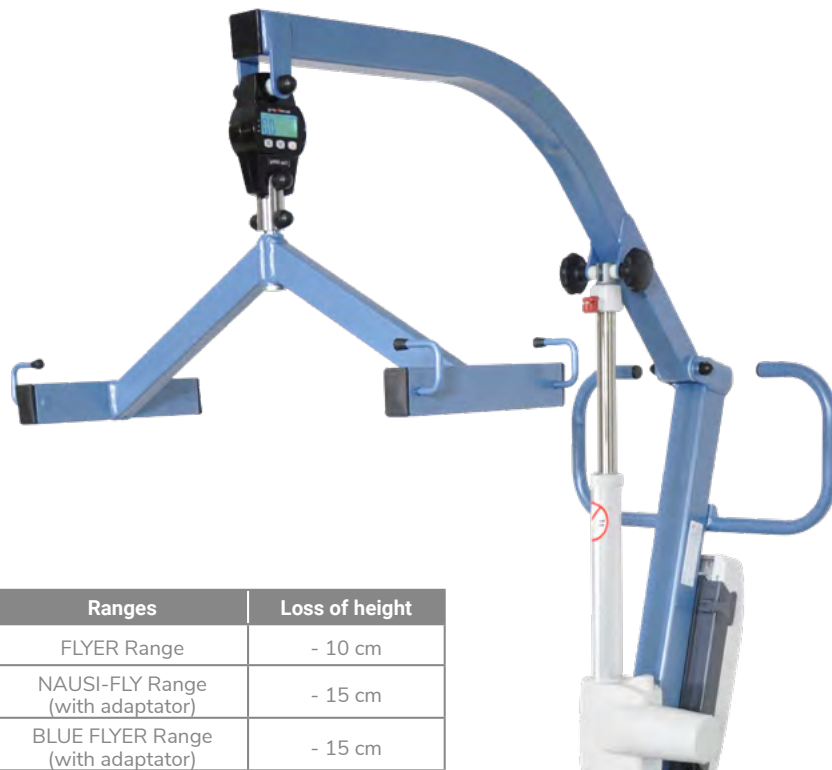

Ranges	Loss of height
FLYER Range	- 10 cm
NAUSI-FLY Range (with adaptator)	- 15 cm
BLUE FLYER Range (with adaptator)	- 15 cm



- 01** TARE: turns the system to zero, even with a weight suspended to the system
- 02** ZERO: used to define the absolute zero
- 03** HOLD: displays a steady weight even if the system oscillates
- 04** UNIT: selects the desired unit of measure

REFERENCES	Weighing System	PESEE-NA-PF
	Adaptator	PESEE-NA-KIT2