

BLUE FLYER

mechanical opening legs



SCAN ME
to learn more

It provides a great simplicity of use to perform secure and comfortable transfers for patients (180 kg weight capacity with an 8 000 N actuator) with its mechanical opening base legs (up to 107 cm) and its proportions which allow a great ease of movement while remaining easy to handle.

WATCH

the BLUE FLYER range product video
on our YOUTUBE channel



YouTube
SCAN
to watch



180 kg



40 kg



620/1800 mm



113 mm




700/1070 mm



360-degree rotation of the spreader bar for a simplified positioning



8 000 N actuator offering an excellent compromise between maximum capacity (180 kg) and lifting speed

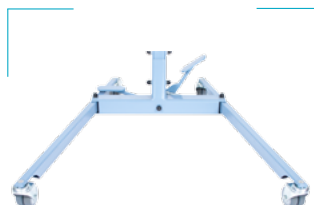

Mechanical opening legs for better access



Reduced height of the front castors to pass under the lowest beds



OPTIONS



Lithium Battery



Weighing System

REFERENCES

with standard battery: FLYER-BLUE-PF
with lithium battery: FLYER-BLUE-LI-PF