



130 kg



WARRANTY
3
YEARS



Decubitus
Level



SCAN ME
to learn more

NAUSIFLOW 2S system

NAUSIFLOW 2S Pump

It operates in alternating/static modes up to 130 kg, is intended for the prevention of very high-risk levels and is used as a care support for stage 3 pressure ulcers.

Its high-speed pump can manage a wide range of mattresses (air height up to 20 cm and width up to 140 cm) and it features a static operating mode to ensure treatment in palliative care.

NAUSIFLOW 100-512 Air Mattress

Its compartmented configuration and its adjusted size (17 x 85 x 200 cm - 13 cm air height - 4 cm foam under-mattress) fits any type of 90 cm width bed (box spring from 85 to 90 cm).


It is equipped with an alternating heel section (4 cells) and a static head section (3 cells) to increase patient comfort.



CPR (cardiopulmonary resuscitation): deflates quickly (less than 20 seconds with a patient) the portion of the mattress which supports the patient's bust



Independent cells in polyurethane on nylon frame

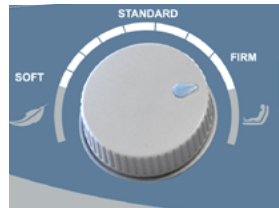

Standard pressure setting: ability to manually adjust the pressure to optimize patient comfort



PALLIATIVE CARE
Palliative care mode: activates the static mode

REFERENCES

NAUSIFLOW 2S PUMP: **NA2S-COMP-PM**
NAUSIFLOW 100-512 MATTRESS: **NA512-MAT85-PM**
NAUSIFLOW 2S SYSTEM: **NA2S512-SYS**