

DEMONSTRATION VIDEO

YouTube
SCAN
to watch



BLUE FLYER PLUS

With Electrical Opening Legs

The second actuator to electrically spread the legs (up to 95 cm) makes it less strenuous for the caregiver who can then devote more attention to the patient.

Equipped with an 8 000 N actuator (180 kg weight capacity), a premium control box and double ball bearing castors for a greater maneuverability, it allows enhanced handling for daily transfers due to its flexibility of use.



KEY FEATURES



360-degree rotation of the spreader bar for a simplified positioning



8 000 N actuator offering a compromise between maximum capacity (180 kg) and lifting speed



Electrical opening legs for ease of use


OPTIONS



A lithium battery (up to 84 cycles - 2 years warranty) is available to enhance the use of the device and a mountable NAUSICAA weighing system to determine the patient's weight.

TECHNICAL DATA

- A. Total length: 119 cm
- B. Minimum lifting height: 60 cm
- C. Maximum lifting height: 181 cm
- D. Total height: 134 cm
- E. Frame height: 11,5 cm
- F. Minimum base width: 66 cm
- G. Maximum base width: 95 cm
- H. Turning diameter: 133 cm

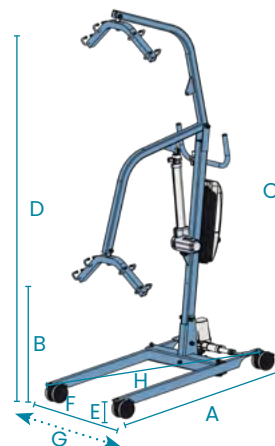

REFERENCES

With standard battery	FLYERBLUE+-E-PF
With lithium battery	FLYERBLUE+-E-LI-PF