


CONVERTISS

Standard Plus



KEY FEATURES

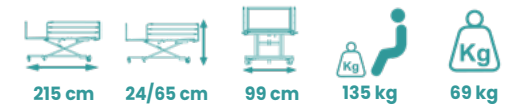

Transport kit



Patient lifting pole in order to easily stand up with a stable point of support



Pair of folding side rails for more safety

PRESENTATION

The CONVERTISS Standard Plus, whose maximum weight capacity is 135 kg, is equipped with 4 actuators allowing 5 different functions: variable height, raising backrest section, adjustable leg rest section, chair position and Trendelenburg position.

Easy to use, it can be split into 4 pieces weighing less than 20 kg and can be assembled without tools.

Its space saving storage facility is the main interest to Nursing, Care Home and Community environments.

The Trendelenburg position (and reverse Trendelenburg) is made to enhance venous return to the heart and prevents pooling of blood in the lower extremities.

Head and foot sections also can be adjusted for the patient's comfort.

Depending on the patient's morphology and pathology of the patient, the height of the bed can be adjusted from 65 cm to 25 cm from the floor to permits an ease of access to the bed.

As an option, in addition to a full-length side rails and a bed lever to enhance the use of the bed, an enlargement kit is available (see page 82) to extend the sleeping width up to 120 cm and increase the length up to 220 cm.

REFERENCE

NAUSIBED-5F